

Method 508 mod
RE-EXTS

Date Extracted 2/4/02
ref

Page 66

Lab #	Description	Sample vol/wt	Drip rate ml/min	Surr add	Spk add	IS add	Final vol.	T.V. pos.	pH	Count																		
<u>LB</u>		<u>500mL</u>		<u>1mL</u>			<u>1mL</u>		<u>7</u>																			
<u>Ref</u>					<u>0.5mL</u>																							
<u>31586</u>								<u>3A</u>																				
<u>31640</u>								<u>5A</u>																				
<u>30646</u>								<u>3B</u>																				
<u>00831</u>								<u>4A</u>																				
<u>00845</u>	<u>2</u>							<u>4B</u>																				
<u>31252</u>								<u>3B</u>																				
<u>31129</u>								<u>2A</u>																				
<u>31132</u>								<u>1B</u>																				
Remarks Reviewed by, <u>CB 5/24/02</u> drip rate must be between 5-15 ml/min for liq-liq extractors drip rate - Extraction type - <u>Septan</u> Solv.lot# <u>V41470</u> NaCl lot# <u>V39618</u> Na2SO4 lot# <u>V39H00</u> <table border="1"> <tr> <td>STANDARD</td> <td>INT STD LOG#</td> </tr> <tr> <td>CALIBRATION</td> <td><u>11013002</u></td> </tr> <tr> <td>REF/SPIKE</td> <td><u>91759</u></td> </tr> <tr> <td>MDL/MDRF</td> <td></td> </tr> <tr> <td>Q. CHECK</td> <td></td> </tr> <tr> <td>INT. STD.</td> <td><u>lot# 61110320</u></td> </tr> <tr> <td>SURROGATE</td> <td><u>020102</u></td> </tr> <tr> <td>LPC</td> <td></td> </tr> <tr> <td>ENDRIN/DDT</td> <td></td> </tr> </table>											STANDARD	INT STD LOG#	CALIBRATION	<u>11013002</u>	REF/SPIKE	<u>91759</u>	MDL/MDRF		Q. CHECK		INT. STD.	<u>lot# 61110320</u>	SURROGATE	<u>020102</u>	LPC		ENDRIN/DDT	
STANDARD	INT STD LOG#																											
CALIBRATION	<u>11013002</u>																											
REF/SPIKE	<u>91759</u>																											
MDL/MDRF																												
Q. CHECK																												
INT. STD.	<u>lot# 61110320</u>																											
SURROGATE	<u>020102</u>																											
LPC																												
ENDRIN/DDT																												

Analyst LF
revision 4/12/01

Verified by my

Reviewed by DN
Date 2/12/02

Sequence Name: C:\HPCHEM\1\SEQUENCE\FEB0802.S

Comment: INST 209; 30Mx0.25umX0.50

Operator: 209 GALOS

Data Path: C:\HPCHEM\1\data\feb0802\

Pre-Seq Cmd:

Post-Seq Cmd:

Method Sections To Run

On A Barcode Mismatch

(X) Full Method

(X) Inject Anyway

() Reprocessing Only

() Don't Inject

Line	Type	Vial	DataFile	Method	Sample Name
------	------	------	----------	--------	-------------

1	Sample	1	DFTPP	GCMS0103	
2	Sample	2	SURR	GCMS0208	
3	Sample	3	CALCHK	GCMS0208	
4	Sample	4	CAL40	GCMS0208	
5	Sample	5	LB	GCMS0208	GC/MS SCAN 2/8
6	Sample	6	REF	GCMS0208	GC/MS SCAN 2/8
7	Sample	7	31586	GCMS0208	GC/MS SCAN 2/4
8	Sample	8	31640	GCMS0208	GC/MS SCAN 2/4
9	Sample	9	30646	GCMS0208	GC/MS SCAN 2/4
10	Sample	10	00831	GCMS0208	GC/MS SCAN 2/4
11	Sample	11	00845	GCMS0208	GC/MS SCAN 2/4
12	Sample	12	31252	GCMS0208	GC/MS SCAN 2/4
13	Sample	13	31129	GCMS0208	GC/MS SCAN 2/4
14	Sample	14	31132	GCMS0208	GC/MS SCAN 2/4
15	Sample	15	LBA	GCMS0208	GC/MS SCAN 1/28
16	Sample	16	REFA	GCMS0208	GC/MS SCAN 1/28
17	Sample	17	1937	GCMS0208	GC/MS SCAN 1/28
18	Sample	18	1938	GCMS0208	GC/MS SCAN 1/28
19	Sample	19	1939	GCMS0208	GC/MS SCAN 1/28
20	Sample	20	1940	GCMS0208	GC/MS SCAN 1/28
21	Sample	21	1941	GCMS0208	GC/MS SCAN 1/28
22	Sample	22	1942	GCMS0208	GC/MS SCAN 1/28
23	Sample	23	1943	GCMS0208	GC/MS SCAN 1/28
24	Sample	24	28049	GCMS0208	GC/MS SCAN 1/28 FB
25	Sample	25	31341	GCMS0208	RE-EXTRACT 1/28
26	Sample	26	DFTPP1	GCMS0208	
27	Sample	27	CAL40A	GCMS0208	
28	Sample	28	LBB	GCMS0208	GC/MS SCAN 2/7
29	Sample	29	REFB	GCMS0208	GC/MS SCAN 2/7
30	Sample	30	2966	GCMS0208	GC/MS SCAN 2/7
31	Sample	31	2967	GCMS0208	GC/MS SCAN 2/7
32	Sample	32	2968	GCMS0208	GC/MS SCAN 2/7
33	Sample	33	2969	GCMS0208	GC/MS SCAN 2/7
34	Sample	34	2970	GCMS0208	GC/MS SCAN 2/7
35	Sample	35	2971	GCMS0208	GC/MS SCAN 2/7
36	Sample	36	2972	GCMS0208	GC/MS SCAN 2/7
37	Sample	37	2973	GCMS0208	GC/MS SCAN 2/7FB
38	Sample	38	3014	GCMS0208	GC/MS SCAN 2/7
39	Sample	39	3015	GCMS0208	GC/MS SCAN 2/7
40	Sample	40	31344	GCMS0208	RE-EXTRACT 2/7
41	Sample	41	31136	GCMS0208	RE-EXTRACT 2/7 FB
42	Sample	42	LBC	GCMS0208	GC/MS SCAN 2/8


Last Modified: Fri Feb 08 16:33:57 2002

Page: 1

DFTPP

Data File : C:\HPCHEM\1\DATA\FEB0802\DFTPP.D
 Acq On : 8 Feb 2002 9:14 am
 Sample :
 Misc :
 MS Integration Params: GOOFY.P
 Method : C:\HPCHEM\1\METHODS\GCMS0208.M (RTE Integrator)
 Title : 209 625/TCLP calibration table

Vial: 1
 Operator: 209 GALOS
 Inst : GC/MS 209
 Multiplr: 1.00

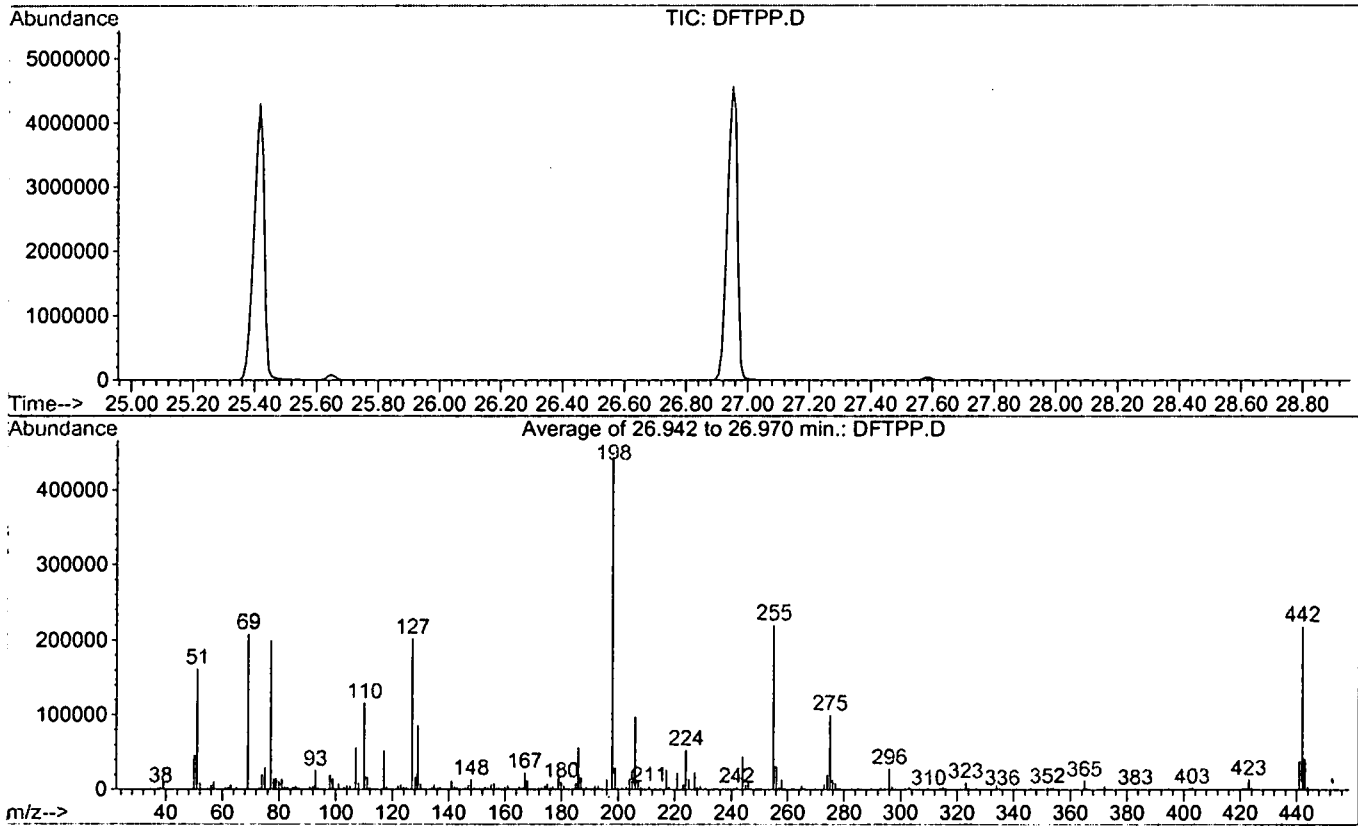

Spectrum Information: Average of 26.942 to 26.970 min.

Target Mass	Rel. to Mass	Lower Limit%	Upper Limit%	Rel. Abn%	Raw Abn	Result Pass/Fail
51	198	30	60	36.4	161730	PASS
68	69	0.00	2	0.0	0	PASS
69	198	0.00	100	46.9	208040	PASS
70	69	0.00	2	0.0	0	PASS
127	198	40	60	45.6	202334	PASS
197	198	0.00	1	0.0	0	PASS
198	198	100	100	100.0	444032	PASS
199	198	5	9	6.7	29625	PASS
275	198	10	30	22.6	100160	PASS
365	198	1	100	2.8	12222	PASS
441	443	0.01	100	92.9	38718	PASS
442	198	40	110	49.3	219072	PASS
443	442	17	23	19.0	41660	PASS