

Method 508 mod

Date Extracted 5-17-02

Page 131

Ref

\$GCMS-18699 / 25728/02

Lab #	Description	Sample vol/wt	Drip rate ml/min	Surr add	Spk add	IS add	Final vol.	T.V. pos.	pH	Count
<u>LB</u>		<u>500ul</u>		<u>1ml</u>			<u>1ml</u>	<u>4B</u>	<u>7</u>	
<u>Rep 1</u>					<u>1ml</u>			<u>4A</u>		
<u>Rep 2</u>					<u>↓</u>			<u>3B</u>		
<u>11424</u>	<u>FB</u>							<u>9B</u>		
<u>11425</u>								<u>2A</u>		
<u>11426</u>								<u>5A</u>		
<u>11427</u>	<u>FB</u>							<u>10A</u>		
<u>11480</u>								<u>9A</u>		
<u>11481</u>								<u>3A</u>		
<u>11482</u>	<u>FB</u>							<u>2B</u>		
<u>11566</u>								<u>10B</u>		
<u>10519</u>	<u>Re-Ext</u>	<u>↓</u>		<u>↓</u>			<u>↓</u>	<u>1A</u>	<u>↓</u>	
Remarks Reviewed by <u>CB 8/23/02</u>										
drip rate must be between 5-15 ml/min for liq-liq extractors drip rate - Extraction type <u>sep fm</u> Solv. lot# <u>X07E07</u> NaCl lot# <u>V39614</u> Na2SO4 lot# <u>V46648</u>										
STANDARD	INT STD LOG#									
CALIBRATION										
REF/SPIKE	<u>2297E</u>									
MDL/MDRF										
Q. CHECK										
INT. STD.										
SURROGATE	<u>TCMX</u>	<u>2284C</u>	<u>5/17/02</u>							
LPC										
ENDRIN/DDT										

Analyst EH
revision 4/12/01

Verified by _____

Reviewed by ON 5/29/02
Date

Method 508 GC/MS scan
Instrument 214

File Directory 020521

Storage Directory _____

Lab #	Filename	Description	Vial #		Inj.	Standards used - Lot #			Remarks	Initials
			front	rear		Calibration	Reference	Surrogate		
	4777P	diff	1							MN
	052101	cov	2			2287i	2287i/2286i			
	057261	lab black	3					2284c	2286i	
	↓ 2	↓ ad	4							
	057741	ref	5			2277c				
	↓ 2	↓ Kb	6							
	↓ 3	↓ ad	7							
	↓ 4	↓ ad	8							
	10579	sample	9							
	11424		10							
	11425		11							
	11426		12							
	11427		13							
	11480		14							
	11481		15							
	11482		16							
	11566		17							
	10812		18							
	11154		19							
	11562		20							
	11568		21							
	11569		22							
	11570		23							
	11571		24							
	11572		25							

Date 05/24/02

Extraction date 05/17/02

Page 11

Sequence Name: C:\MSDchem\1\sequence\020521.S
Comment: GCMS scan
Operator: Nefliu
Data Path: C:\MSDchem\1\DATA\020521\
Pre-Seq Cmd:
Post-Seq Cmd:


Method Sections To Run On A Barcode Mismatch
(X) Full Method (X) Inject Anyway
() Reprocessing Only () Don't Inject

Line	Type	Vial	DataFile	Method	Sample Name
1	Sample	1	DFTPP	508FAST	
2	Sample	2	0521C1	508SCAN1	ccv
3	Sample	3	0517LB1	508SCAN1	05/17 lab blk kb
4	Sample	4	0517LB2	508SCAN1	05/17 lab blk ad
5	Sample	5	0517REF1	508SCAN1	052/17 ref kb
6	Sample	6	0517REF2	508SCAN1	052/17 ref kb
7	Sample	7	0517REF3	508SCAN1	052/17 ref ad
8	Sample	8	0517REF4	508SCAN1	052/17 ref ad
9	Sample	9	10519	508SCAN1	sample
10	Sample	10	11424	508SCAN1	sample
11	Sample	11	11425	508SCAN1	sample
12	Sample	12	11426	508SCAN1	sample
13	Sample	13	11427	508SCAN1	sample
14	Sample	14	11480	508SCAN1	sample
15	Sample	15	11481	508SCAN1	sample
16	Sample	16	11482	508SCAN1	sample
17	Sample	17	11566	508SCAN1	sample
18	Sample	18	10812	508SCAN1	sample
19	Sample	19	11154	508SCAN1	sample
20	Sample	20	11567	508SCAN1	sample
21	Sample	21	11568	508SCAN1	sample
22	Sample	22	11569	508SCAN1	sample
23	Sample	23	11570	508SCAN1	sample
24	Sample	24	11571	508SCAN1	sample
25	Sample	25	11572	508SCAN1	sample

DFTPP

Data File : C:\MSDCHEM\1\DATA\020521\DFTPP.D
 Acq On : 21 May 2002 8:20 am
 Sample :
 Misc :
 MS Integration Params: rteint.e
 Method : C:\MSDCHEM\1\METHODS\SCAN020521.M (RTE Integrator)
 Title : Method 508 GC/MS Scan

Vial: 1
 Operator: Nefliu
 Inst : Instrume
 Multiplr: 1.00

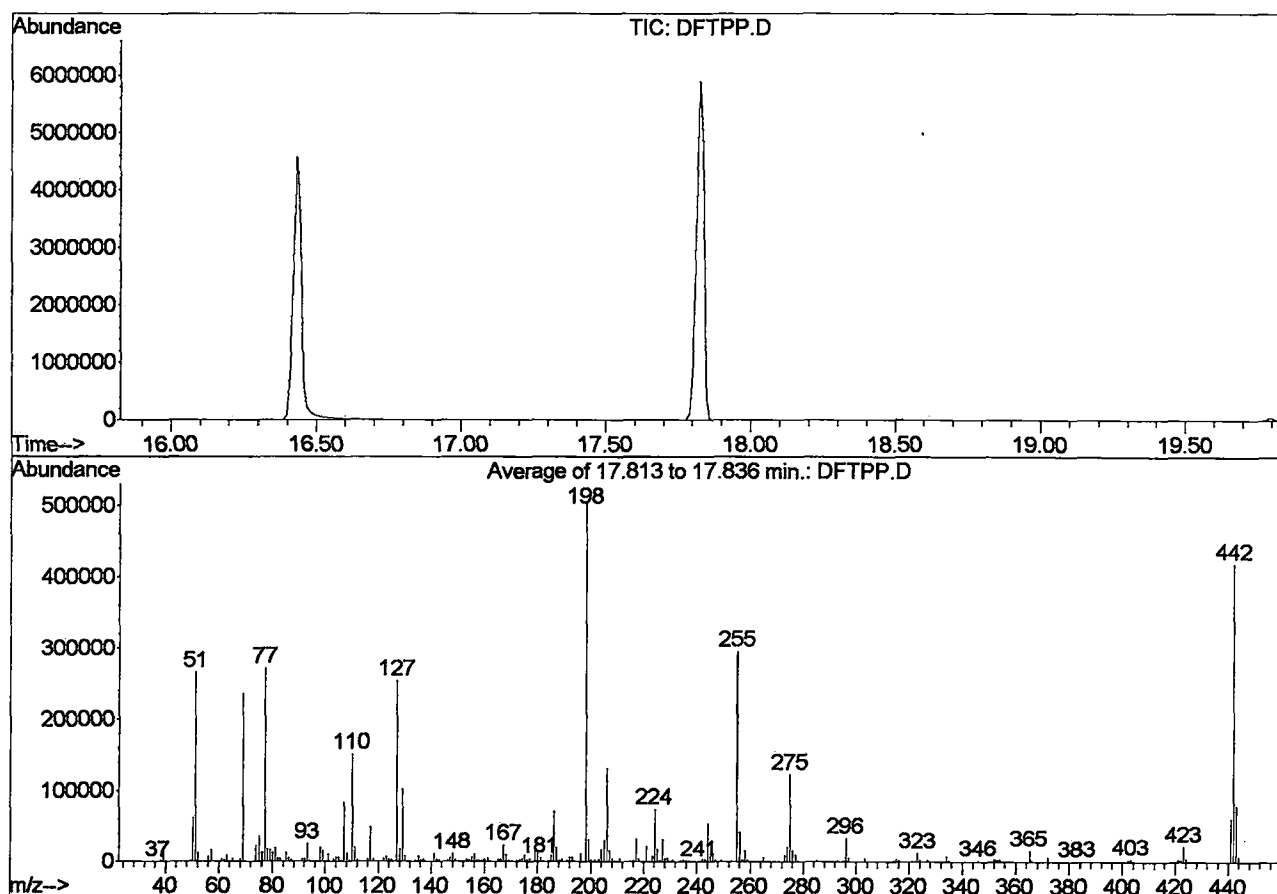

Spectrum Information: Average of 17.813 to 17.836 min.

Target Mass	Rel. to Mass	Lower Limit%	Upper Limit%	Rel. Abn%	Raw Abn	Result Pass/Fail
51	198	30	60	52.7	266368	PASS
68	69	0.00	2	1.6	3881	PASS
69	198	0.00	100	46.7	236416	PASS
70	69	0.00	2	0.5	1263	PASS
127	198	40	60	50.3	254464	PASS
197	198	0.00	1	0.0	0	PASS
198	198	100	100	100.0	505749	PASS
199	198	5	9	6.4	32432	PASS
275	198	10	30	24.5	123808	PASS
365	198	1	99	3.3	16669	PASS
441	443	1	100	78.4	62418	PASS
442	198	40	110	82.8	418816	PASS
443	442	17	23	19.0	79653	PASS