

Method 508 MOD.

Date Extracted 4/26/02

Page 84

Ref

*with another
re-extended*

Lab #	Description	Sample vol/wt	Drip rate ml/min	Surr add	Spk add	IS add	Final vol.	T.V. pos.	pH	Count		
L.B.		500 ml		1 ml			1 ml	9B	7			
Ref 1					1 ml			7B				
Ref 2					1 ml			7A				
9151	Re-EXT.	↓						3B				
9154	Re-EXT.	415 ml						3A				
9155	Re-EXT.	490 ml						4B				
9331	Re-EXT.	500 ml						8B				
9332	Re-EXT.							8A				
10062								4A				
10063								10B				
10064								9A				
10065								10A				
Remarks <i>Original batch - many low surrogate recoveries</i>												
drip rate must be between 5-15 ml/min for liq-liq extractors drip rate - Extraction type - <i>SEP. Fun.</i> Solv. lot# <i>X05E10</i> NaCl lot# <i>V 39614</i> Na2SO4 lot# <i>V 46648</i>												
STANDARD	INT STD LOG#											
CALIBRATION												
REF/SPIKE	<i>2279E</i>											
MDL/MDRF												
Q. CHECK												
INT. STD.	<i>042602</i>											
SURROGATE	<i>GCMs</i>											
LPC												
ENDRIN/DDT												

Analyst A.D.
Vision 4/12/01

Verified by _____

Reviewed by _____
Date

Sequence Name: C:\HPCHEM\1\SEQUENCE\MAY0302.S

Comment:

Operator:

Data Path: C:\HPCHEM\1\DATA\MAY0302\

Pre-Seq Cmd:

Post-Seq Cmd:

Method Sections To Run

(X) Full Method

() Reprocessing Only

On A Barcode Mismatch

(X) Inject Anyway

() Don't Inject

Line	Type	Vial	DataFile	Method	Sample Name
1	Sample	1	DFTPP	GCMS0501	
2	Sample	2	CAL20	GCMS0501	
3	Sample	3	LB	GCMS0501	GC/MS SCAN 4/24
4	Sample	4	REF1	GCMS0501	GC/MS SCAN 4/24
5	Sample	5	REF2	GCMS0501	GC/MS SCAN 4/24
6	Sample	6	8837	GCMS0501	RE-EXT.
7	Sample	7	8838	GCMS0501	RE-EXT.
8	Sample	8	8934	GCMS0501	RE-EXT.
9	Sample	9	9579	GCMS0501	GC/MS SCAN 4/24
10	Sample	10	9580	GCMS0501	GC/MS SCAN 4/24
11	Sample	11	9581	GCMS0501	GC/MS SCAN 4/24 FB
12	Sample	12	9585	GCMS0501	GC/MS SCAN 4/24
13	Sample	13	9586	GCMS0501	GC/MS SCAN 4/24
14	Sample	14	9587	GCMS0501	GC/MS SCAN 4/24 FB
15	Sample	15	LBA	GCMS0501	GC/MS SCAN 4/22
16	Sample	16	REFA1	GCMS0501	GC/MS SCAN 4/22
17	Sample	17	REFA2	GCMS0501	GC/MS SCAN 4/22
18	Sample	18	9517	GCMS0501	GC/MS SCAN 4/22
19	Sample	19	9518	GCMS0501	GC/MS SCAN 4/22
20	Sample	20	9519	GCMS0501	GC/MS SCAN 4/22
21	Sample	21	9520	GCMS0501	GC/MS SCAN 4/22
22	Sample	22	8271	GCMS0501	RE-EXT. FB
23	Sample	23	8430	GCMS0501	RE-EXT.
24	Sample	24	8485	GCMS0501	RE-EXT.
25	Sample	25	8487	GCMS0501	RE-EXT.
26	Sample	26	9446	GCMS0501	RE-EXT.
27	Sample	1	DFTPP2	GCMS0501	
28	Sample	100	CAL40	GCMS0501	
29	Sample	27	LBB	GCMS0501	GC/MS SCAN 4/26
30	Sample	28	REFB1	GCMS0501	GC/MS SCAN 4/26
31	Sample	29	REFB2	GCMS0501	GC/MS SCAN 4/26
32	Sample	30	9151	GCMS0501	RE-EXT.
33	Sample	31	9154	GCMS0501	RE-EXT.
34	Sample	32	9155	GCMS0501	RE-EXT.
35	Sample	33	9331	GCMS0501	RE-EXT.
36	Sample	34	9332	GCMS0501	RE-EXT.
37	Sample	35	10062	GCMS0501	GC/MS SCAN 4/26
38	Sample	36	10063	GCMS0501	GC/MS SCAN 4/26
39	Sample	37	10064	GCMS0501	GC/MS SCAN 4/26
40	Sample	38	10065	GCMS0501	GC/MS SCAN 4/26
41	Sample	39	LBC	GCMS0501	GC/MS SCAN 4/23
42	Sample	40	REFC1	GCMS0501	GC/MS SCAN 4/23
43	Sample	41	REFC2	GCMS0501	GC/MS SCAN 4/23


Sequence Name: C:\HPCHEM\1\SEQUENCE\MAY0302.S

Line	Type	Vial	DataFile	Method	Sample Name
44	Sample	42 9568	GCMS0501	GC/MS	SCAN 4/23 FB
45	Sample	43 9569	GCMS0501	GC/MS	SCAN 4/23
46	Sample	44 9570	GCMS0501	GC/MS	SCAN 4/23
47	Sample	45 9571	GCMS0501	GC/MS	SCAN 4/23
48	Sample	46 9572	GCMS0501	GC/MS	SCAN 4/23 FB
49	Sample	47 9574	GCMS0501	GC/MS	SCAN 4/23
50	Sample	48 9575	GCMS0501	GC/MS	SCAN 4/23
51	Sample	49 9576	GCMS0501	GC/MS	SCAN 4/23
52	Sample	50 9577	GCMS0501	GC/MS	SCAN 4/23
53	Sample	51 9578	GCMS0501	GC/MS	SCAN 4/23

DFTPP

Data File : C:\HPCHEM\1\DATA\MAY0302\DFTPP2.D
 Acq On : 4 May 2002 3:34 pm
 Sample :
 Misc :
 MS Integration Params: GOOFY.P
 Method : C:\HPCHEM\1\METHODS\GCMS0501.M (RTE Integrator)
 Title : 209 625/TCLP calibration table

Vial: 1
 Operator:
 Inst : 209
 Multiplr: 1.00

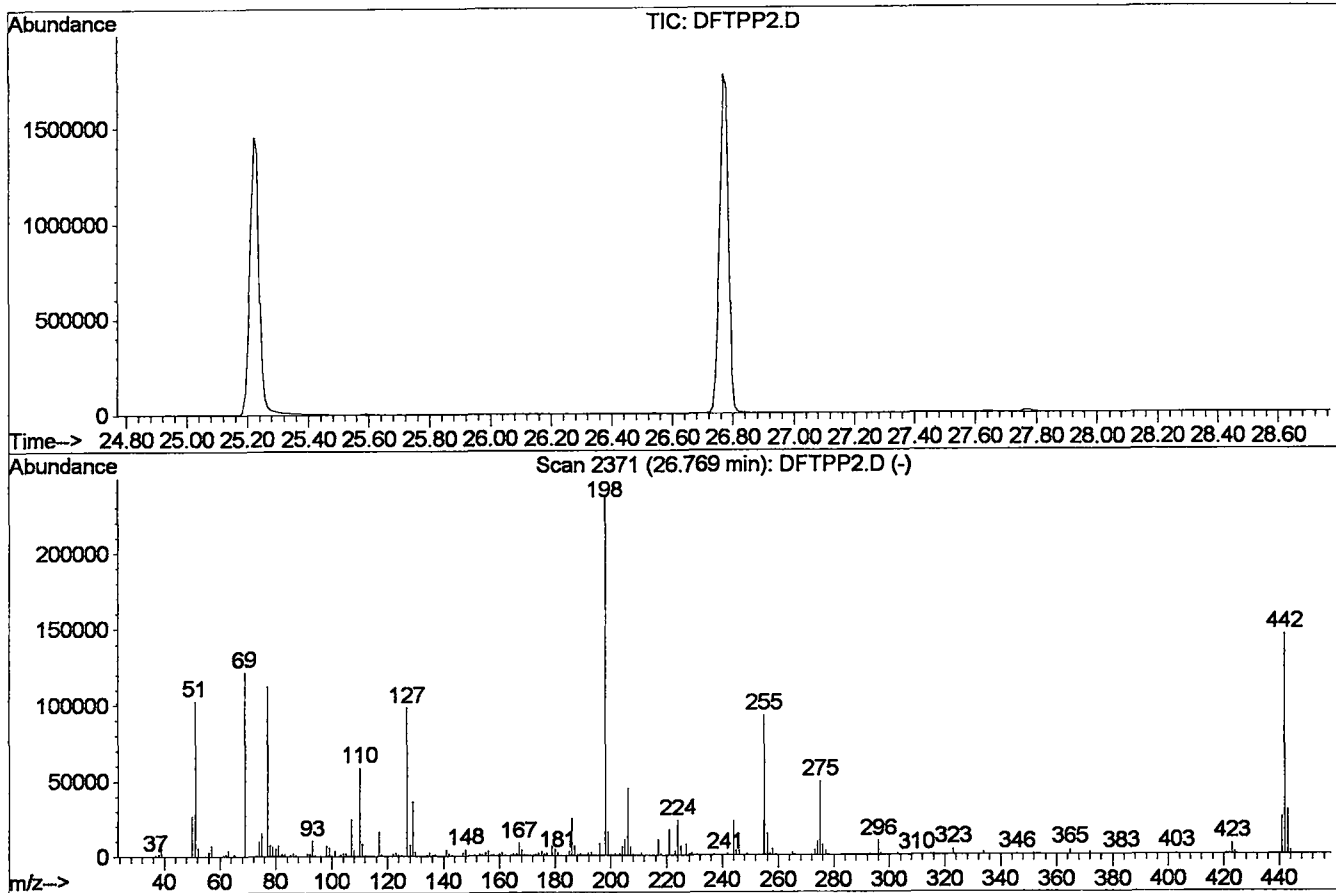

Spectrum Information: Scan 2371

Target Mass	Rel. to Mass	Lower Limit%	Upper Limit%	Rel. Abn%	Raw Abn	Result Pass/Fail
51	198	30	60	43.3	102688	PASS
68	69	0.00	2	0.0	0	PASS
69	198	0.00	100	51.4	121912	PASS
70	69	0.00	2	0.0	0	PASS
127	198	40	60	41.4	98256	PASS
197	198	0.00	1	0.0	0	PASS
198	198	100	100	100.0	237056	PASS
199	198	5	9	6.5	15505	PASS
275	198	10	30	20.7	49064	PASS
365	198	1	100	1.4	3390	PASS
441	443	0.01	100	84.8	24704	PASS
442	198	40	110	61.1	144896	PASS
443	442	17	23	20.1	29120	PASS