

. . . 508 modified

Date Extracted 11/23/01

Page 14

46CMS-15714

Lab #	Description	Sample vol/wt	Drip rate ml/min	Surr add	Spk add	IS add	Final vol.	T.V. pos.	pH	Count
LB	✓	500mL		1mL			1mL		7	
REF	✓				0.5mL					
28600	✓									
28595	✓									
28596	✓									
28473	-HD ✓									
28479	-HD ✓									
28477	-HD ✓									
28598	✓									
28475	✓ HD ✓									
28597	✓									
28599	✓									
28652	✓									
28594	✓									
28593	✓									
28649	✓									
28651	✓									
28650	✓ broken ml resistant -11/28									
28665	✓									
28666	✓									
28667	✓									
28668	✓									
Reviewed by, <u>QA 5/11/02</u>										
Remarks: 28649, 28650, 28651: after buffer: pH=2 → + ~3M NaOH 8N.										
drip rate must be between 5-15 ml/min for liq-liq extractors										
drip rate -										
Extraction type - <u>sepm</u>		Solv. lot#		NaCl lot# <u>V11549</u>		Na2SO4 lot# <u>V32603</u>				
STANDARD	INT STD LOG#									
CALIBRATION										
REF/SPIKE	<u>91759</u>									
MDL/MDRF										
Q. CHECK										
INT. STD.										
SURROGATE	<u>24000</u>									
LPC										
ENDRIN/DDT										

Analyst LF/MN
revision 4/12/01

Verified by _____

Reviewed by DW
Date 12/7/11

Sequence Name: C:\HPCHEM\1\SEQUENCE\NOV2901.S

Comment: INST 209; 30Mx0.25umX0.50

Operator: 209 GALOS

Data Path: C:\HPCHEM\1\data\nov2901\

Pre-Seq Cmd:

Post-Seq Cmd:

Method Sections To Run

On A Barcode Mismatch

☒ Full Method☒ Inject Anyway☐ Reprocessing Only☐ Don't Inject

Line Type	Vial	DataFile	Method	Sample Name
1 Sample	1	DFTPP	GCMS1126	
2 Sample	2	CALCHK	GCMS1126	25
3 Sample	24	28468	GCMS1126	11/21 scan
4 Sample	25	28431	GCMS1126	11/21 scan
5 Sample	26	28434	GCMS1126	11/21 scan
6 Sample	27	28427	GCMS1126	11/21 scan
7 Sample	28	28428	GCMS1126	11/21 scan
8 Sample	31	28469	GCMS1126	11/21 scan
9 Sample	3	CALI25	GCMS1126	calib 25
10 Sample	32	28275	GCMS1126	11/21 scan
11 Sample	33	28263	GCMS1126	11/21 scan
12 Sample	34	28591	GCMS1126	11/23 scan
13 Sample	35	28666	GCMS1126	11/23 scan
14 Sample	36	28593	GCMS1126	11/23 scan
15 Sample	37	28595	GCMS1126	11/23 scan
16 Sample	38	28475	GCMS1126	11/23 scan
17 Sample	39	28600	GCMS1126	11/23 scan
18 Sample	40	28596	GCMS1126	11/23 scan
19 Sample	41	28652	GCMS1126	11/23 scan
20 Sample	42	28665	GCMS1126	11/23 scan
21 Sample	43	28599	GCMS1126	11/23 scan
22 Sample	44	28668	GCMS1126	11/23 scan
23 Sample	45	28667	GCMS1126	11/23 scan
24 Sample	46	28597	GCMS1126	11/23 scan
25 Sample	47	28598	GCMS1126	11/23 scan
26 Sample	48	28649	GCMS1126	11/23 scan
27 Sample	1	DFTPP1	GCMS1126	
28 Sample	3	CALI25A	GCMS1126	calib 25
29 Sample	32	28275A	GCMS1126	11/21 scan
30 Sample	49	28651	GCMS1126	11/23 scan
31 Sample	51	LBD	NP112601	lb 11/2 gc/ms
32 Sample	52	LBE	NP112601	lb 11/2 gc/ms
33 Sample	53	26079	NP112601	11/2 gc/ms
34 Sample	54	26707	NP112601	11/2 gc/ms
35 Sample	55	26208	NP112601	11/2 gc/ms
36 Sample	56	26206	NP112601	11/2 gc/ms
37 Sample	57	26251	NP112601	11/2 gc/ms
38 Sample	58	26080	NP112601	11/2 gc/ms
39 Sample	59	26133	NP112601	11/2 gc/ms
40 Sample	60	26163	NP112601	11/2 gc/ms
41 Sample	61	26164	NP112601	11/2 gc/ms
42 Sample	62	26165	NP112601	11/2 gc/ms
43 Sample	63	26127	NP112601	11/2 gc/ms


Last Modified: Fri Nov 30 10:30:40 2001

Page: 1

DFTPP

Data File : C:\HPCHEM\1\DATA\NOV2901\DFTPP.D
 Acq On : 29 Nov 2001 8:42 am
 Sample :
 Misc :
 MS Integration Params: RTEINT.P
 Method : C:\HPCHEM\1\METHODS\GCMS1126.M (RTE Integrator)
 Title : 209 625/TCLP calibration table

Vial: 1
 Operator: 209 GALOS
 Inst : GC/MS 209
 Multiplr: 1.00

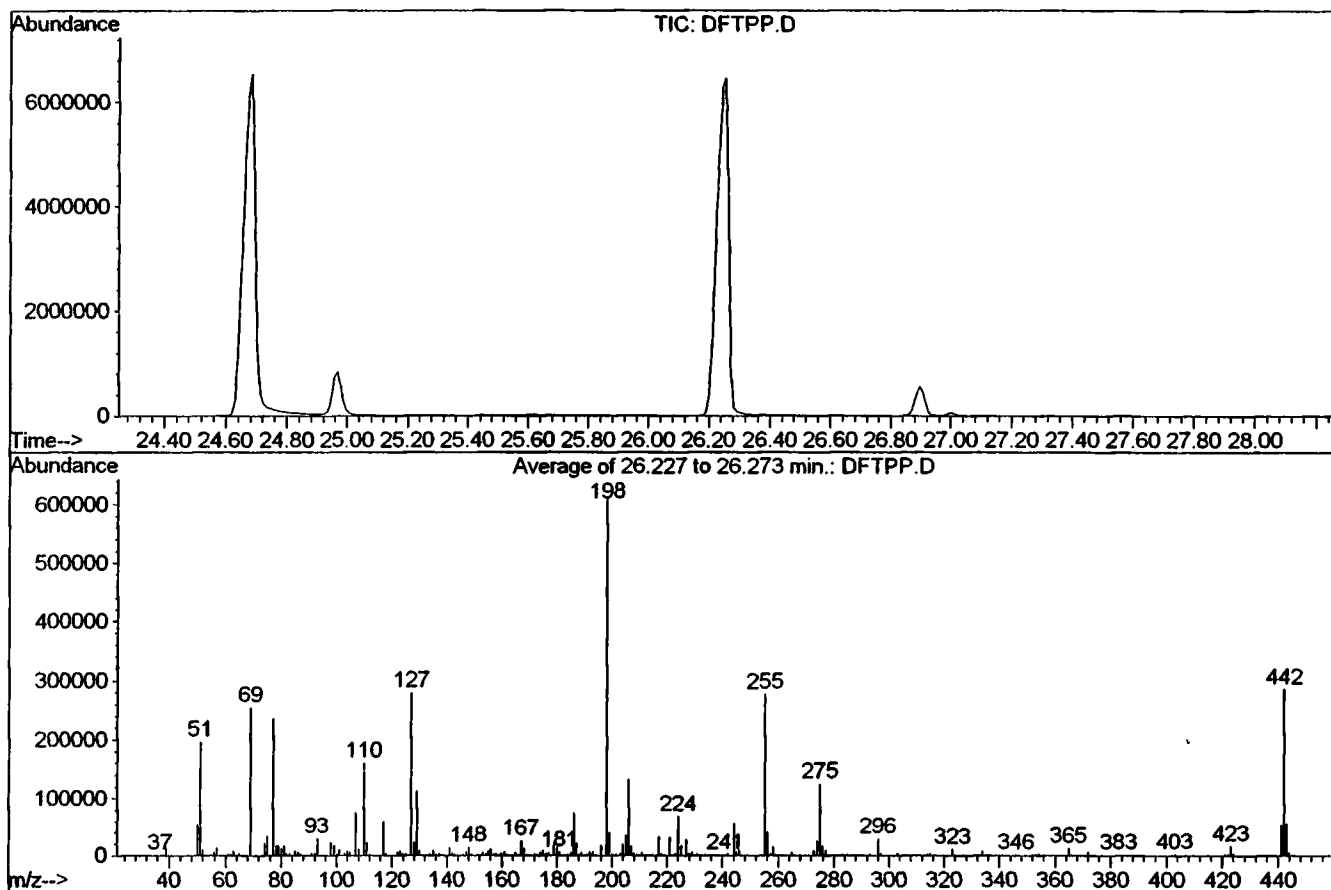

Spectrum Information: Average of 26.227 to 26.273 min.

Target Mass	Rel. to Mass	Lower Limit%	Upper Limit%	Rel. Abn%	Raw Abn	Result Pass/Fail
51	198	30	60	32.1	196867	PASS
68	69	0.00	2	0.0	0	PASS
69	198	0.00	100	41.7	255692	PASS
70	69	0.00	2	0.0	29	PASS
127	198	40	60	45.7	280255	PASS
197	198	0.00	1	0.0	0	PASS
198	198	100	100	100.0	613224	PASS
199	198	5	9	6.6	40720	PASS
275	198	10	30	20.3	124768	PASS
365	198	1	100	2.0	12568	PASS
441	443	0.01	100	96.0	53009	PASS
442	198	40	110	46.7	286411	PASS
443	442	17	23	19.3	55196	PASS