

Method 508 modified

Date Extracted 12/5/01

Page 41

GCMS14132

REF

Lab #	Description	Sample vol/wt	Drip rate ml/min	Surr add	Spk add	IS add	Final vol.	T.V. pos.	pH	Count		
✓ LB ✓		500 µL		1 µL			1 µL		7			
✓ REF ✓	1 Bottle (P) Bottle				0.5 µL							
✓ 29457 ✓												
✓ 29458 ✓												
✓ 29290 ✓												
✓ 29289 ✓												
✓ 29291 ✓												
✓ 29459 ✓												
✓ 29460 ✓												
✓ 29461 ✓												
✓ 29462 ✓												
✓ 29298 ✓												
✓ 29297 ✓												
✓ 29345 ✓												
✓ 29294 ✓	(P) Bottle											
	not needed											
	biochem sample											

Remarks

All the samples were Acidic including the Ref & LB because buffer was made with the wrong solute (monobasic dipotassium crystal). REF + LB were re-extracted using Dibasic potassium powder and all the samples were given 2mL 6M NaOH to bring the pH to 7.

drip rate must be between 5-15 ml/min for liq-liq extractors

drip rate -

Extraction type - Sup fun

Solv. lot#

NaCl lot#

V11599

Na2SO4 lot#

V39H00

STANDARD	INT STD LOG#											
CALIBRATION												
REF/SPIKE	<u>91759</u>											
MDL/MDRF												
Q. CHECK												
INT. STD.												
SURROGATE	<u>GCMS 112601</u>											
LPC												
ENDRIN/DDT												

Reviewed by,

4/5 3/14/02

Analyst Heath
revision 4/12/01

Verified by _____

Reviewed by
Date 1/4/02


Sequence Name: C:\HPCHEM\1\SEQUENCE\DCE1001.S

Line Type	Vial	DataFile	Method	Sample Name
44 Sample	5	29456	GCMS1126	gc/ms scan 12/4
45 Sample	6	29455	GCMS1126	gc/ms scan 12/4
46 Sample	7	29295	GCMS1126	gc/ms scan 12/4
47 Sample	8	29235	GCMS1126	gc/ms scan 12/4
48 Sample	9	29296	GCMS1126	gc/ms scan 12/4
49 Sample	10	29294	GCMS1126	gc/ms scan 12/4
50 Sample	11	29291	GCMS1126	gc/ms scan 12/4
51 Sample	12	29229	GCMS1126	gc/ms scan 12/4
52 Sample	13	29285	GCMS1126	gc/ms scan 12/4
53 Sample	14	29293	GCMS1126	gc/ms scan 12/4
54 Sample	15	LBE	GCMS1126	gc/ms scan 12/5
55 Sample	16	REFE	GCMS1126	gc/ms scan 12/5
56 Sample	17	29457	GCMS1126	gc/ms scan 12/5
57 Sample	18	29461	GCMS1126	gc/ms scan 12/5
58 Sample	19	29460	GCMS1126	gc/ms scan 12/5
59 Sample	20	29594	GCMS1126	gc/ms scan 12/5
60 Sample	21	29345	GCMS1126	gc/ms scan 12/5
61 Sample	22	29462	GCMS1126	gc/ms scan 12/5
62 Sample	23	29297	GCMS1126	gc/ms scan 12/5
63 Sample	24	29289	GCMS1126	gc/ms scan 12/5
64 Sample	25	29298	GCMS1126	gc/ms scan 12/5

DFTPP

Data File : C:\HPCHEM\1\DATA\DEC1201\DFTPP.D
 Acq On : 12 Dec 2001 10:47 am
 Sample :
 Misc :
 MS Integration Params: GOOFY.P
 Method : C:\HPCHEM\1\METHODS\GCMS1126.M (RTE Integrator)
 Title : 209 625/TCLP calibration table

Vial: 1
 Operator: 209 GALOS
 Inst : GC/MS 209
 Multiplr: 1.00

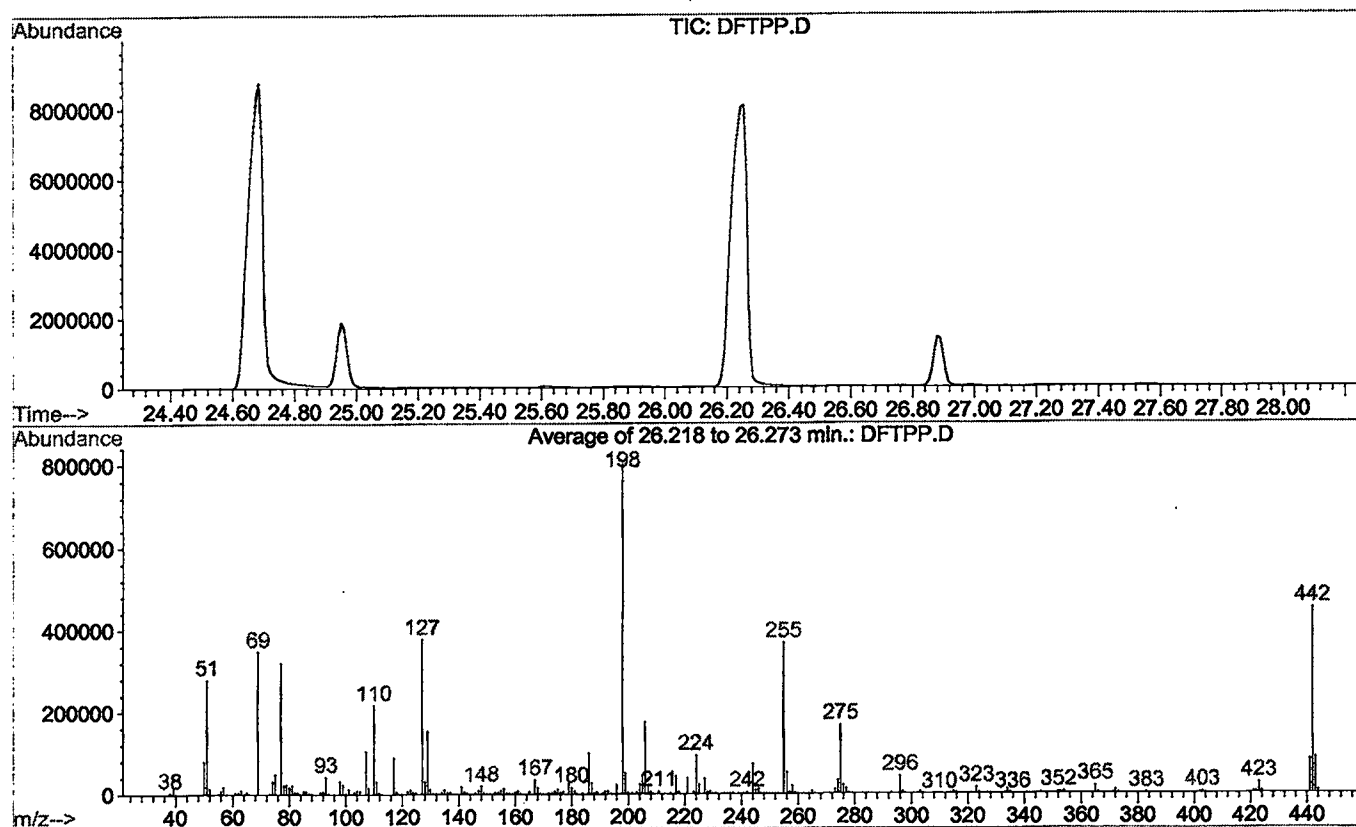

Spectrum Information: Average of 26.218 to 26.273 min.

Target Mass	Rel. to Mass	Lower Limit%	Upper Limit%	Rel. Abn%	Raw Abn	Result Pass/Fail
51	198	30	60	35.1	280730	PASS
68	69	0.00	2	0.0	0	PASS
69	198	0.00	100	43.5	348389	PASS
70	69	0.00	2	0.0	77	PASS
127	198	40	60	47.2	377585	PASS
197	198	0.00	1	0.0	0	PASS
198	198	100	100	100.0	800229	PASS
199	198	5	9	6.6	52431	PASS
275	198	10	30	21.0	167777	PASS
365	198	1	100	2.4	19596	PASS
441	443	0.01	100	94.6	82423	PASS
442	198	40	110	56.3	450514	PASS
443	442	17	23	19.3	87131	PASS