

Method 508 MODDate Extracted 5-2-02Page Batch
#18554

REF

Lab #	Description	Sample vol/wt	Drip rate ml/min	Surr add	Spt add	IS add	Final vol.	T.V. pos.	pH	Count		
	LAB BLANK	500 ml		1 ml			1 ml	2A	7			
	REF 1				1 ml			2B				
	REF 2				1 ml			3A				
	10431							3B				
	10432							4A				
	10433							4B				
	10434							6A				
	10435							6B				
	10436							7A				
	10437							7B				
Remarks												
Reviewed by, <u>S/B 8/22/02</u>												
drip rate must be between 5-15 ml/min for liq-liq extractors												
drip rate -												
Extraction type - <u>SEP. FUNNEL</u> Solv. lot# <u>X05E10</u> NaCl lot# <u>V 39614</u> Na2SO4 lot# <u>V 46648</u>												
STANDARD	INT STD LOG#											
CALIBRATION	<u>HPD 5020</u>											
REF SPIKE	<u>2279E 4/19/02</u>											
MDL/MDRF												
Q. CHECK												
INT. STD.	<u>056702 05/002</u>											
SURROGATE	<u>IST 440-2284 5/4/02</u>											
LPC												
ENDRIN/DDT												

Analyst B.N
revision 4/12/01Verified by myReviewed by DN
Date 5/20/02

Sequence Name: C:\HPCHEM\1\SEQUENCE\MAY1002.S

Comment:

Operator:

Data Path: C:\HPCHEM\1\DATA\MAY1002\

Pre-Seq Cmd:

Post-Seq Cmd:

Method Sections To Run

On A Barcode Mismatch

(X) Full Method

(X) Inject Anyway

() Reprocessing Only

() Don't Inject

Line	Type	Vial	DataFile	Method	Sample Name
1	Sample	1	DFTPP	GCMS0507	
2	Sample	2	CAL40	GCMS0507	
3	Sample	3	LB	GCMS0510	GC/MS SCAN 5/9
4	Sample	4	REF1	GCMS0510	GC/MS SCAN 5/9
5	Sample	5	REF2	GCMS0510	GC/MS SCAN 5/9
6	Sample	6	9398	GCMS0510	GC/MS SCAN 5/9
7	Sample	7	10605	GCMS0510	GC/MS SCAN 5/9
8	Sample	8	9683	GCMS0510	GC/MS SCAN 5/9
9	Sample	9	9764	GCMS0510	GC/MS SCAN 5/9
10	Sample	10	9765	GCMS0510	GC/MS SCAN 5/9
11	Sample	11	9766	GCMS0510	GC/MS SCAN 5/9
12	Sample	12	9770	GCMS0510	GC/MS SCAN 5/9
13	Sample	13	LBA	GCMS0507	GC/MS SCAN 5/6
14	Sample	14	REFA1	GCMS0507	GC/MS SCAN 5/6
15	Sample	15	REFA2	GCMS0507	GC/MS SCAN 5/6
16	Sample	16	10772	GCMS0507	GC/MS SCAN 5/6
17	Sample	17	10773	GCMS0507	GC/MS SCAN 5/6
18	Sample	18	10774	GCMS0507	GC/MS SCAN 5/6
19	Sample	19	10775	GCMS0507	GC/MS SCAN 5/6 fb
20	Sample	20	10776	GCMS0507	GC/MS SCAN 5/6
21	Sample	21	10777	GCMS0507	GC/MS SCAN 5/6
22	Sample	22	10778	GCMS0507	GC/MS SCAN 5/6
23	Sample	23	10779	GCMS0507	GC/MS SCAN 5/6 fb
24	Sample	1	DFTPP1	GCMS0507	
25	Sample	24	CAL80	GCMS0507	
26	Sample	25	LBB	GCMS0507	gc/ms scan 5/2
27	Sample	26	REFB1	GCMS0507	gc/ms scan 5/2
28	Sample	27	REFB2	GCMS0507	gc/ms scan 5/2
29	Sample	28	10431	GCMS0507	gc/ms scan 5/2
30	Sample	29	10432	GCMS0507	gc/ms scan 5/2
31	Sample	30	10433	GCMS0507	gc/ms scan 5/2
32	Sample	31	10434	GCMS0507	gc/ms scan 5/2
33	Sample	32	10435	GCMS0507	gc/ms scan 5/2
34	Sample	33	10436	GCMS0507	gc/ms scan 5/2
35	Sample	34	10437	GCMS0507	gc/ms scan 5/2
36	Sample	35	LBC	GCMS0507	gc/ms scan 5/2
37	Sample	36	REFC1	GCMS0507	gc/ms scan 5/2
38	Sample	37	REFC2	GCMS0507	gc/ms scan 5/2
39	Sample	38	10515	GCMS0507	gc/ms scan 5/2
40	Sample	39	10516	GCMS0507	gc/ms scan 5/2
41	Sample	40	10517	GCMS0507	gc/ms scan 5/2
42	Sample	41	10518	GCMS0507	gc/ms scan 5/2 fb
43	Sample	42	10519	GCMS0507	gc/ms scan 5/2


Sequence Name: C:\HPCHEM\1\SEQUENCE\MAY1002.S

Line Type	Vial	DataFile	Method	Sample Name
44 Sample	43 10520	GCMS0507	gc/ms	scan 5/2
45 Sample	44 10521	GCMS0507	gc/ms	scan 5/2
46 Sample	45 10522	GCMS0507	gc/ms	scan 5/2 fb

DFTPP

Data File : C:\HPCHEM\1\DATA\MAY1002\DFTPP1.D
 Acq On : 11 May 2002 11:32 am
 Sample :
 Misc :
 MS Integration Params: GOOFY.P
 Method : C:\HPCHEM\1\METHODS\GCMS0510.M (RTE Integrator)
 Title : 209 625/TCLP calibration table

Vial: 1
 Operator:
 Inst : 209
 Multiplr: 1.00

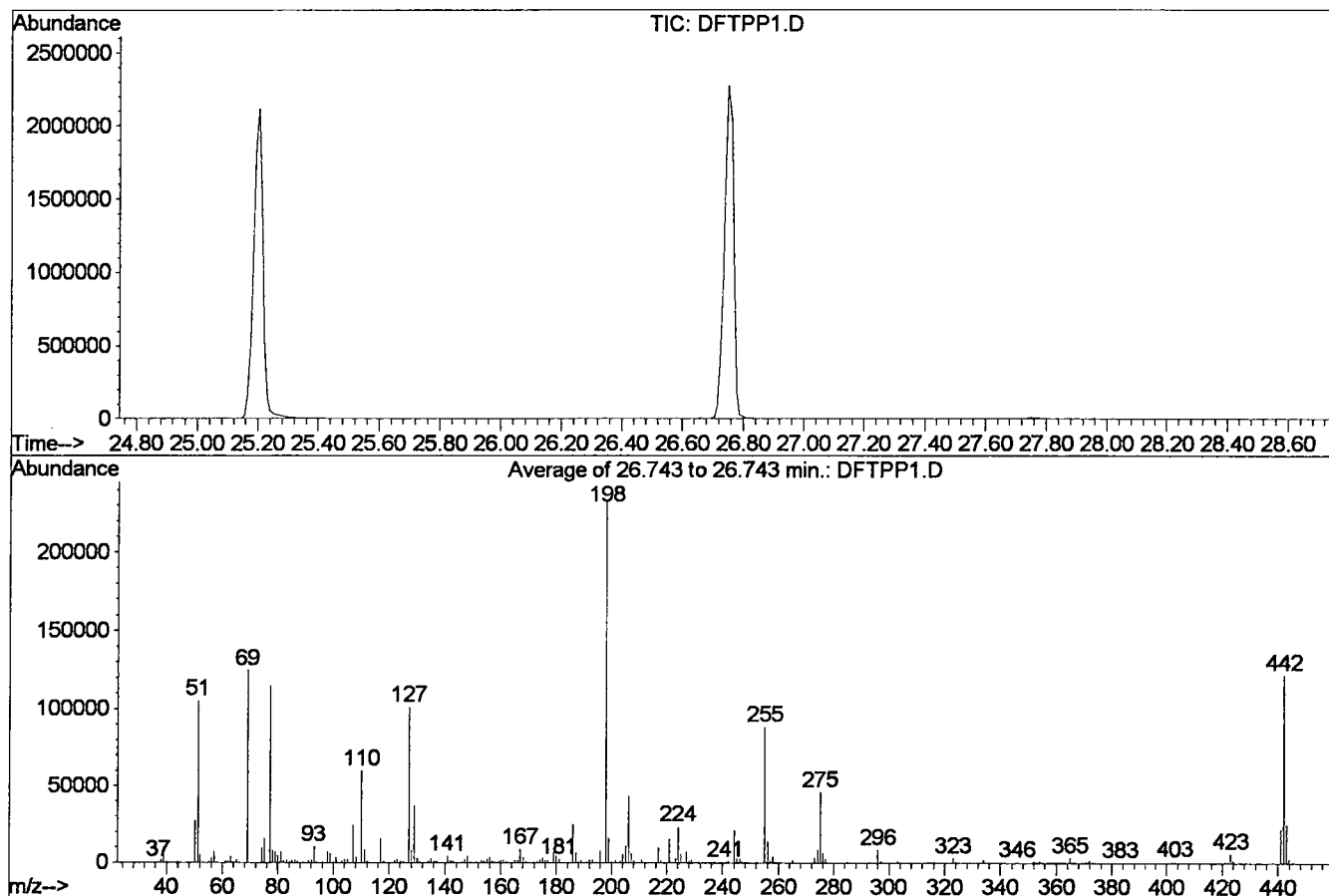

Spectrum Information: Average of 26.743 to 26.743 min.

Target Mass	Rel. to Mass	Lower Limit%	Upper Limit%	Rel. Abn%	Raw Abn	Result Pass/Fail
51	198	30	60	45.1	105304	PASS
68	69	0.00	2	0.0	0	PASS
69	198	0.00	100	53.4	124680	PASS
70	69	0.00	2	0.0	0	PASS
127	198	40	60	43.4	101296	PASS
197	198	0.00	1	0.0	0	PASS
198	198	100	100	100.0	233600	PASS
199	198	5	9	6.8	15950	PASS
275	198	10	30	19.7	46104	PASS
365	198	1	100	1.4	3203	PASS
441	443	0.01	100	87.3	21480	PASS
442	198	40	110	52.3	122144	PASS
443	442	17	23	20.1	24592	PASS