

Method 508 mod.

Date Extracted 3-29-02

Page 35

Lab #	Description	Sample vol/wt	Drip rate ml/min	Surr add	Spk add	IS add	Final vol.	T.V. pos.	pH	Count
103		500ml		1ml	1ml		1ml	6A	7	
7416								2A		
7417								2B		
7418								1A		
7419								4A		
7420								4B		
7421								1B		
7422								7A		
7423	FB	↓	↓				↓	7B	↓	
								8A		
								8B		
Remarks										
Reviewed by, <u>CB 8/14/02</u>										
drip rate must be between 5-15 ml/min for liq-liq extractors drip rate - Extraction type <u>SEP Fur</u> Spk lot# <u>XDIE19</u> NaCl lot# <u>V39618</u> Na2SO4 lot# <u>V46623</u>										
STANDARD	INT STD LOG#									
CALIBRATION										
REF SPIKE	<u>2193</u>									
MOD/MORF										
Q CHECK										
INT. STD.										
SURROGATE										
LPC										
ENDRN/DOF										

Analyst EA
Revision 4/12/01

Verified by _____

Reviewed by DN
Date 7/26/02

508SCAN

Method 525.2

Instrument 213

File Directory 041902

Storage Directory

Lab #	Filename	Description	Vial #		Inj.	Standards used - Lot #			Int. Std.	Remarks	Initials
			front	rear		ul	Calibration	Reference			
	DFTA0419	DFTPP	1		1	21106					TB
	C80A0419	C80	2			020412					
	SURA0419	Surrogate	3							DDD	
	SURTO419	1 g	4							TCMX, 3/25/02	
	LBA0419	LB	5					2286b			
	REFA0419	Ref 1	6					1			
	SURTO419	Surrogate	25							TCMX, 4/9/02	
	REFB0419	Ref 2	7					2286b			
	MCA0419	MCA0419	8								
	7416	Sample	9					2286b			
	7417		10							based on	
	7418		11							cal of	
	7419		12							4/12/02	
	7420		13								
	7421		14								
	7422		15								
	7423		16								
	MCB0419	MC	17								
	SURX419	Surrogate	18							TCMX, 3/25/02	
	SURY419	1 g	19							TCMX, 4/9/02	
	C80B0419	C80	20			020412					
File DFTA0419 was accidentally erased. DFTPP was reviewed prior to erasure and it passed. TB											

Date 4/19/02

Extraction date 3/29/02

Page

Sequence Name: C:\MSDCHEM\1\SEQUENCE\041902.S

Comment:

Operator: Bohlk

Data Path: C:\MSDCHEM\1\DATA\041902\

Pre-Seq Cmd:

Post-Seq Cmd:


Method Sections To Run On A Barcode Mismatch

(X) Full Method (X) Inject Anyway

() Reprocessing Only () Don't Inject

Line Type	Vial	DataFile	Method	Sample Name
1 Sample	1	DFTA0419	508SCAN	DFTPP
2 Sample	2	C80A0419	508SCAN	C80
3 Sample	3	SURA0419	508SCAN	SURROGATE, DDD
4 Sample	4	SURT0419	508SCAN	SURROGATE, TCMX 3/25/02
5 Sample	5	LBA0419	508SCAN	LB 3/29/02
6 Sample	6	REFA0419	508SCAN	REF1 3/29/02
7 Sample	25	SURJ0419	508SCAN	SURROGATE, TCMX 4/9/02
8 Sample	7	REFB0419	508SCAN	REF2 3/29/02
9 Sample	8	MCA0419	508SCAN	MC
10 Sample	9	7416	508SCAN	SAMPLE 7416 3/29/02
11 Sample	10	7417	508SCAN	SAMPLE 7416 3/29/02
12 Sample	11	7418	508SCAN	SAMPLE 7418 3/29/02
13 Sample	12	7419	508SCAN	SAMPLE 7419 3/29/02
14 Sample	13	7420	508SCAN	SAMPLE 7420 3/29/02
15 Sample	14	7421	508SCAN	SAMPLE 7421 3/29/02
16 Sample	15	7422	508SCAN	SAMPLE 7422 3/29/02
17 Sample	16	7423	508SCAN	SAMPLE 7423 3/29/02
18 Sample	17	MCB0419	508SCAN	MC
19 Sample	18	SURX419	508SCAN	SURR, TCMX 3/25/02
20 Sample	19	SURY419	508SCAN	SURR, TCMX 4/19/02
21 Sample	20	C80B0419	508SCAN	C80

File : C:\MSDCHEM\1\DATA\041902\DFTA0419.D
 Operator : Bohlk
 Acquired : 19 Apr 2002 12:18 pm using AcqMethod 508SCAN
 Instrument : Instrumen
 Sample Name: DFTPP
 Misc Info :
 Vial Number: 1

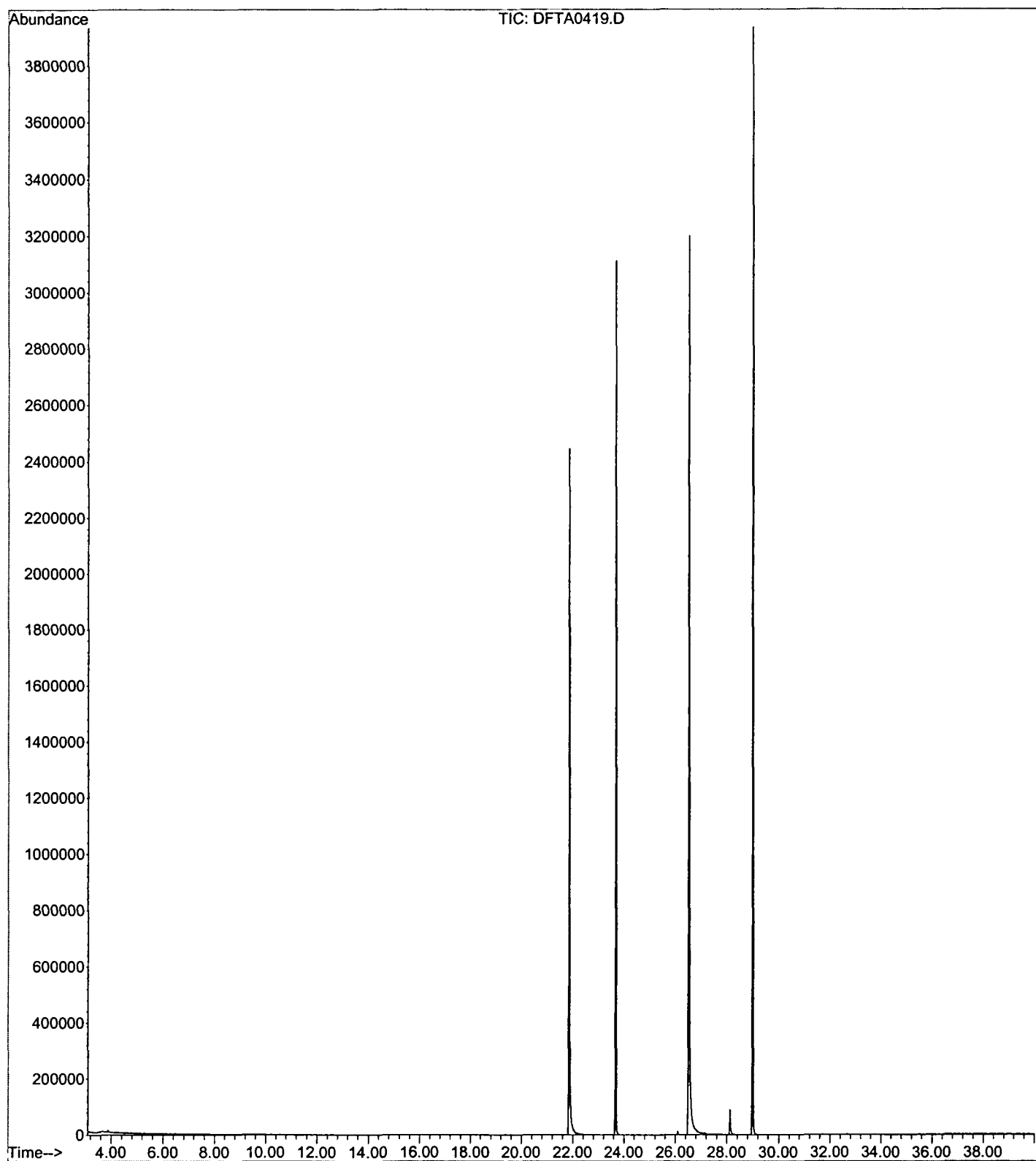

DFTPP

Data File : C:\MSDCHEM\1\DATA\041902\DFTA0419.D
Acq On : 19 Apr 2002 12:18 pm
Sample : DFTPP
Misc :
MS Integration Params: BENZO.P

Vial: 1
Operator: Bohlk
Inst : Instrumen
Multiplr: 1.00

Method : C:\MSDCHEM\1\5973N\5080402BBC.M (RTE Integrator)
Title : Quantitation For Method 508gcscan

Spectrum Information: Scan 3506

Target Mass	Rel. to Mass	Lower Limit%	Upper Limit%	Rel. Abn%	Raw Abn	Result Pass/Fail
51	198	30	60	46.1	173248	PASS
68	69	0.00	2	1.5	2562	PASS
69	198	0.00	100	46.9	176384	PASS
70	69	0.00	2	0.7	1255	PASS
127	198	40	60	51.7	194304	PASS
197	198	0.00	1	0.0	0	PASS
198	198	100	100	100.0	376192	PASS
199	198	5	9	6.8	25752	PASS
275	198	10	30	19.2	72160	PASS
365	198	1	100	2.5	9235	PASS
441	443	0.01	100	90.9	35792	PASS
442	198	40	110	56.9	213952	PASS
443	442	17	23	18.4	39392	PASS

DFTA0419.D 5080402BBC.M Fri Apr 19 13:06:35 2002