

Method 508
Ge/ms

Date Extracted 10/10/01

Page 22

Lab #	Description	Sample vol/wt	Drip rate ml/min	Surr add	Spk add	IS add	Final vol.	T.V. pos.	pH	Count
	Lab 6th	500ml	—	1		—	2M		7	
	Ref 10X									
	MDL 10X									
	PR1 1248									
	24215 ✓									2
	24216 ✓									3
	24217 ✓									4
	24218 ✓									5
	24219 ✓									6
	24220 ✓									7
	24236 ✓									8
	24237 ✓									9
	24238 ✓									10
	24287 ✓									11
	24288 ✓									12
	24289 ✓									13
	24290 ✓									14
	24291 ✓									15
	24292 ✓									16
	24293 ✓									17
Reviewed by <u>95</u> <u>1/24/02</u>										
Remarks										
drip rate must be between 5-15 ml/min for liq-liq extractors										
drip rate -										
Extraction type -		Solv. lot#		NaCl lot#		Na2SO4 lot#				
STANDARD	INT STD LOG#									
CALIBRATION										
REF/SPIKE	10X 1678a Perms 1504									
MDL/MDRF										
Q. CHECK										
INT. STD.										
SURROGATE	2064c									
LPC										
ENDRIN/DDT										

Analyst mg/AD/KM
revision 4/12/01

Verified by _____

Reviewed by DN
Date 11/7/01

Sequence Name: C:\HPCHEM\1\SEQUENCE\OCT3101.S
 Comment: INST 209; 30Mx0.25umX0.50
 Operator: 209 GALOS
 Data Path: C:\HPCHEM\1\data\oct3101\
 Pre-Seq Cmd:
 Post-Seq Cmd:

Method Sections To Run On A Barcode Mismatch
 (X) Full Method (X) Inject Anyway
 () Reprocessing Only () Don't Inject

Line Type	Vial	DataFile	Method	Sample Name
-----------	------	----------	--------	-------------

1 Sample	1	DFTPP	NP101901	
2 Sample	2	LB	NP101901	10/10
3 Sample	3	24236	NP101901	gc/ms unknown
4 Sample	4	24237	NP101901	gc/ms unknown
5 Sample	5	24238	NP101901	gc/ms unknown
6 Sample	6	24290	NP101901	gc/ms unknown
7 Sample	7	24291	NP101901	gc/ms unknown
8 Sample	8	24287	NP101901	gc/ms unknown
9 Sample	9	24288	NP101901	gc/ms unknown
10 Sample	10	24289	NP101901	gc/ms unknown
11 Sample	11	24292	NP101901	gc/ms unknown
12				

reinject
reinject

Sequence Name: C:\HPCHEM\1\SEQUENCE\NOV0101.S
Comment: INST 209; 30Mx0.25umX0.50
Operator: 209 GALOS
Data Path: C:\HPCHEM\1\data\nov0101\
Pre-Seq Cmd:
Post-Seq Cmd:

Method Sections To Run	On A Barcode Mismatch
(X) Full Method	(X) Inject Anyway
() Reprocessing Only	() Don't Inject

Line Type	Vial	DataFile	Method	Sample Name
1 Sample	1	DFTPP	NP101901	
2 Sample	2	LB	NP101901	10/10
3 Sample	3	24236	NP101901	gc/ms unknown
4 Sample	2	24215	NP101901	gc/ms unknown
5 Sample	3	24216	NP101901	gc/ms unknown
6 Sample	4	24217	NP101901	gc/ms unknown
7 Sample	5	24218	NP101901	gc/ms unknown
8 Sample	6	24219	NP101901	gc/ms unknown
9 Sample	7	24220	NP101901	gc/ms unknown
10 Sample	8	LBA	NP101901	10/11
11 Sample	9	24333	NP101901	gc/ms unknown
12 Sample	10	24334	NP101901	gc/ms unknown
13 Sample	11	24335	NP101901	gc/ms unknown
14 Sample	12	24395	NP101901	gc/ms unknown
15 Sample	13	24396	NP101901	gc/ms unknown

DFTPP

Data File : C:\HPCHEM\1\DATA\OCT3101\DFTPP.D

Acq On : 31 Oct 2001 4:07 pm

Sample :

Misc :

MS Integration Params: RTEINT.P

Method : C:\HPCHEM\1\METHODS\NP103001.M (RTE Integrator)


Title : 209 625/TCLP calibration table

Vial: 1

Operator: 209 GALOS

Inst : GC/MS 209

Multiplr: 1.00

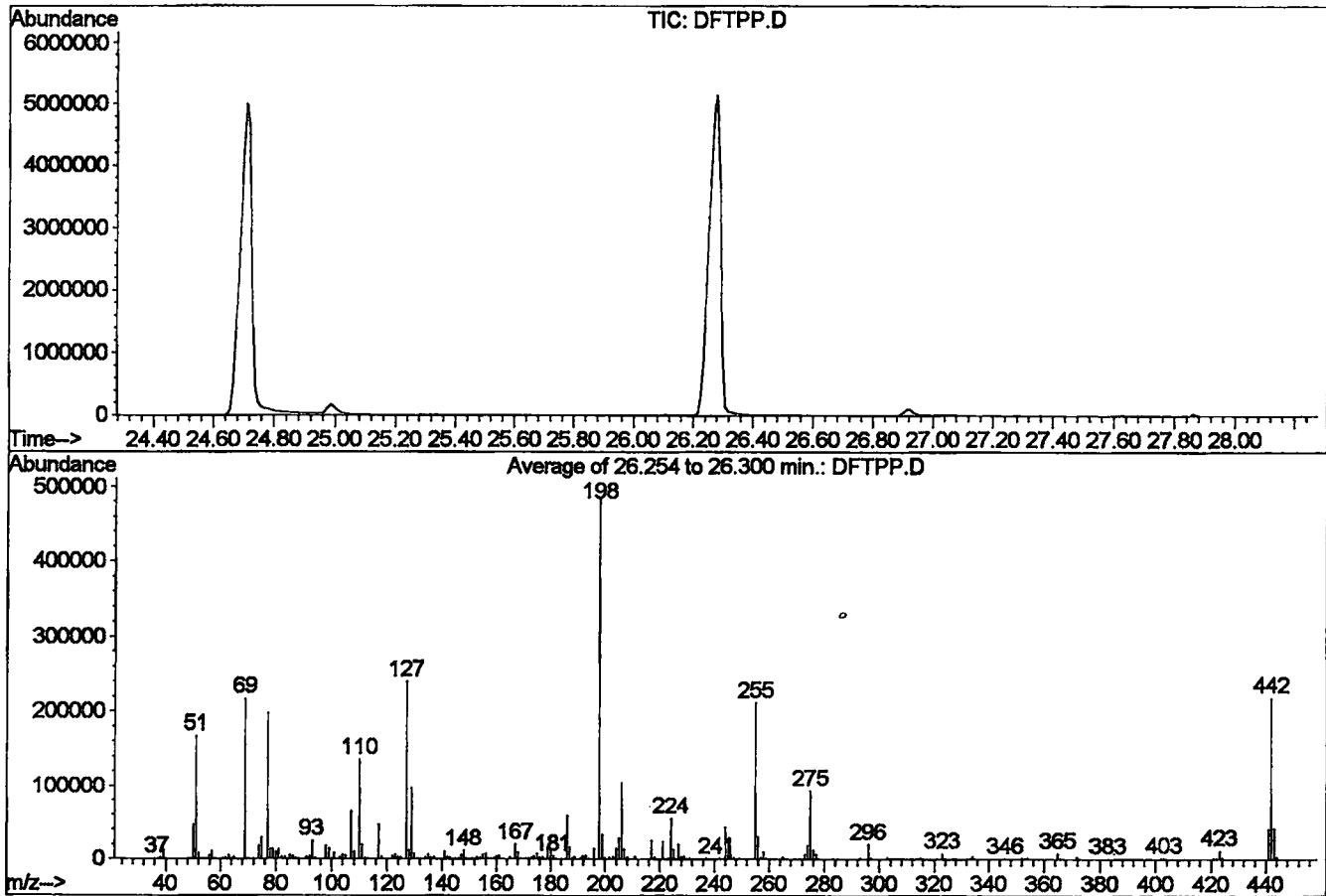

Spectrum Information: Average of 26.254 to 26.300 min.

Target Mass	Rel. to Mass	Lower Limit%	Upper Limit%	Rel. Abn%	Raw Abn	Result Pass/Fail
51	198	30	60	34.4	167296	PASS
68	69	0.00	2	0.0	0	PASS
69	198	0.00	100	44.7	217529	PASS
70	69	0.00	2	0.0	0	PASS
127	198	40	60	49.3	239789	PASS
197	198	0.00	1	0.0	0	PASS
198	198	100	100	100.0	486513	PASS
199	198	5	9	6.8	33113	PASS
275	198	10	30	19.3	94123	PASS
365	198	1	100	1.7	8260	PASS
441	443	0.01	100	98.5	40885	PASS
442	198	40	110	45.0	218923	PASS
443	442	17	23	19.0	41523	PASS

DFTPP.D NP103001.M Thu Nov 01 10:29:34 2001