


File : C:\MSDCHEM\1\DATA\050102\DFTA0501.D
 Operator : Bohlk / Inst 213
 Acquired : 1 May 2002 12:12 pm using AcqMethod 050102
 Instrument : Instrumen
 Sample Name: DFTPP/ENDRIN
 Misc Info :
 Vial Number: 1

Initial Cal

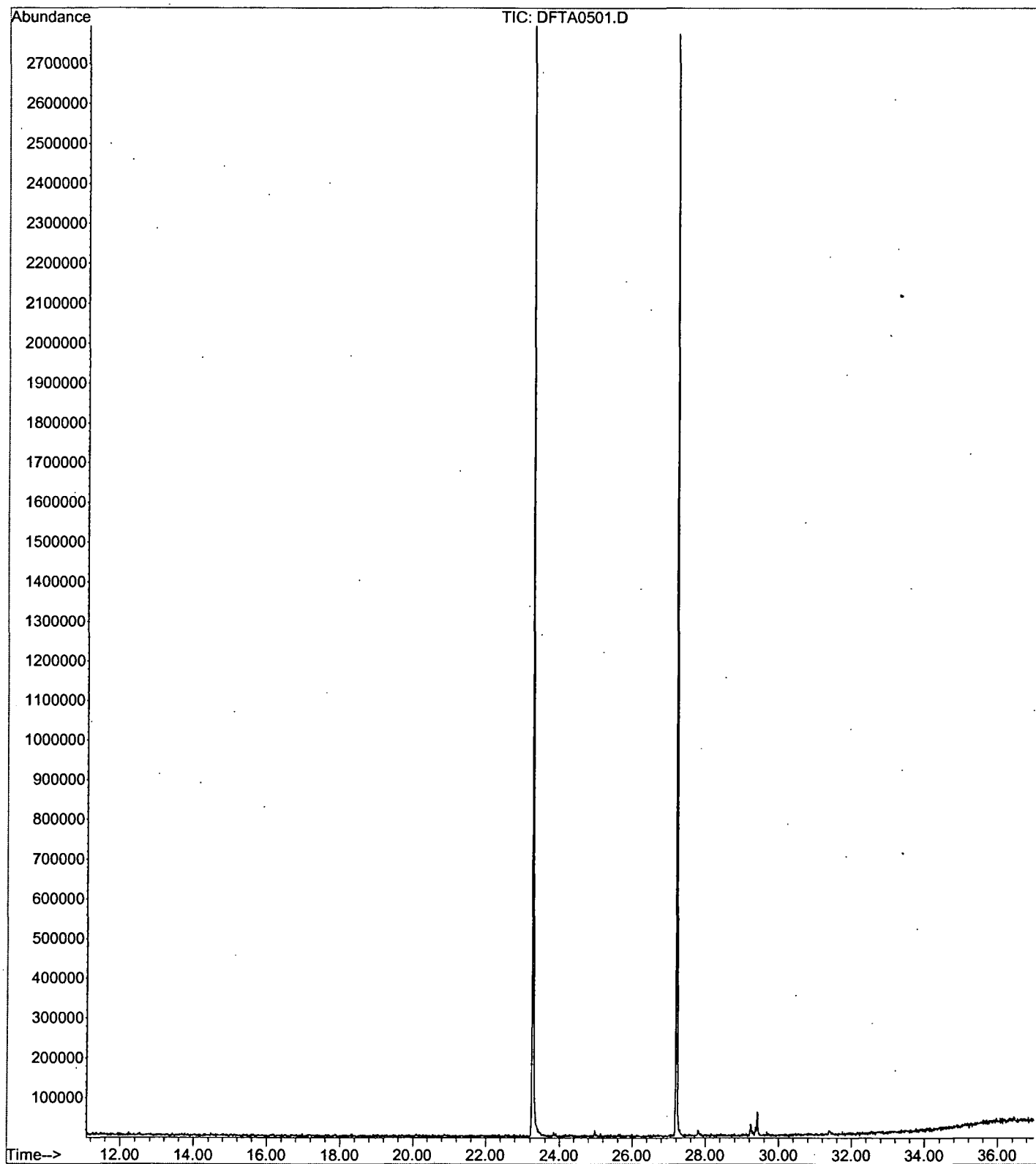

DFTPP

Data File : C:\MSDCHEM\1\DATA\050102\DFTA0501.D

Acq On : 1 May 2002 12:12 pm

Sample : DFTPP/ENDRIN

Misc :

MS Integration Params: rteint.p

Vial: 1

Operator: Bohlk / Inst 213

Inst : Instrumen

Multiplr: 1.00

Method : C:\MSDCHEM\1\5973N\050102.M (RTE Integrator)

Title : Quantitation For Method 525.2

Spectrum Information: Average of 23.185 to 23.439 min.

Target Mass	Rel. to Mass	Lower Limit%	Upper Limit%	Rel. Abn%	Raw Abn	Result Pass/Fail
51	198	10	80 ✓	34.6	15226	PASS
68	69	0.00	2 ✓	1.4	248	PASS
69	198	0.00	100 ✓	39.7	17451	PASS
70	69	0.00	2 ✓	0.4	69	PASS
127	198	10	80 ✓	55.7	24507	PASS
197	198	0.00	2 ✓	0.3	136	PASS
198	198	50	100 ✓	100.0 ✓	43983	PASS
199	198	5	9 ✓	6.4	2796	PASS
275	198	10	60 ✓	18.7	8209	PASS
365	198	1	100 ✓	1.9	827	PASS
441	443	0.01	100 ✓	74.0	3957	PASS
442	198	50	100 ✓	65.3	28737	PASS
443	442	15	24 ✓	18.6	5347	PASS

DFTA0501.D 050102.M

Thu May 02 15:22:41 2002

Area Percent Report

Data File : C:\MSDCHEM\1\DATA\050102\DTA0501.D
Acq On : 1 May 2002 12:12 pm
Sample : DFTPP/ENDRIN
Misc :
MS Integration Params: rteint.p

Vial: 1
Operator: Bohlk / Inst 213
Inst : Instrumen
Multiplr: 1.00

Method : C:\MSDCHEM\1\5973N\050102.M (RTE Integrator)
Title : Quantitation For Method 525.2
Smoothing : ON
Sampling : 1
Start Thrs: 0.2
Stop Thrs : 0
Filtering: 5
Min Area: 5000 Area counts
Max Peaks: 10
Peak Location: TOP

If leading or trailing edge < 100 prefer < Baseline drop else tangent >
Peak separation: 5


Signal : TIC

peak #	R.T. min	first scan	max scan	last scan	PK TY	peak height	corr. area	corr. % max.	% of total
1	27.203	1773	1778	1786	rM	2773114	5574874	100.00%	94.967%
2	27.802	1840	1844	1852	rM 4	14448	42022	0.75%	0.716%
3	29.416	1997	2022	2031	rM 8	60349	253421	4.55%	4.317%

Sum of corrected areas: 5870317

> 5.0 % ✓
Endrin
Breakdown

File : C:\MSDCHEM\1\DATA\050102\DDTA0501.D
Operator : Bohlk / Inst 213
Acquired : 1 May 2002 12:58 pm using AcqMethod 050102
Instrument : Instrumen
Sample Name: DDT
Misc Info :
Vial Number: 2

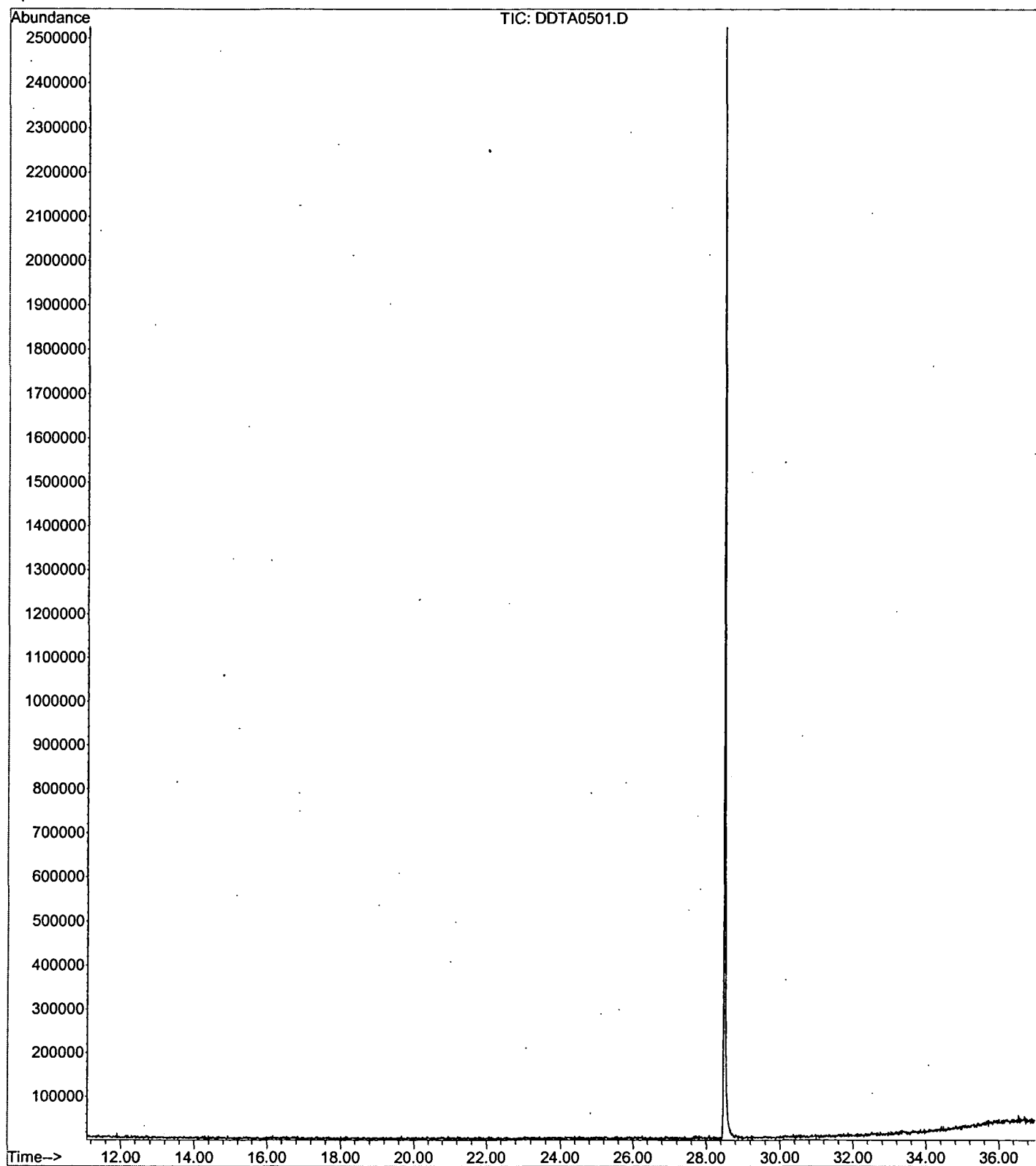

Area Percent Report

J) Data File : C:\MSDCHEM\1\DATA\050102\DDTA0501.D Vial: 2
 Acq On : 1 May 2002 12:58 pm Operator: Bohlk / Inst 213
 Sample : DDT Inst : Instrumen
 Misc : Multiplr: 1.00
 MS Integration Params: rteint.p

Method : C:\MSDCHEM\1\5973N\050102.M (RTE Integrator)
 Title : Quantitation For Method 525.2
 Smoothing : ON Filtering: 5
 Sampling : 1 Min Area: 5000 Area counts
 Start Thrs: 0.2 Max Peaks: 10
 Stop Thrs : 0 Peak Location: TOP

If leading or trailing edge < 100 prefer < Baseline drop else tangent >
 Peak separation: 5

Signal : TIC

peak #	R.T. min	first scan	max scan	last scan	PK TY	peak height	corr. area	corr. % max.	% of total
1	28.491	1914	1920	1948	rM	2519948	5345302	100.00%	100.000%

Sum of corrected areas: 5345302

2% Breakdown